
SYSTEMAT

ORAL IRRIGATOR **EL25921**

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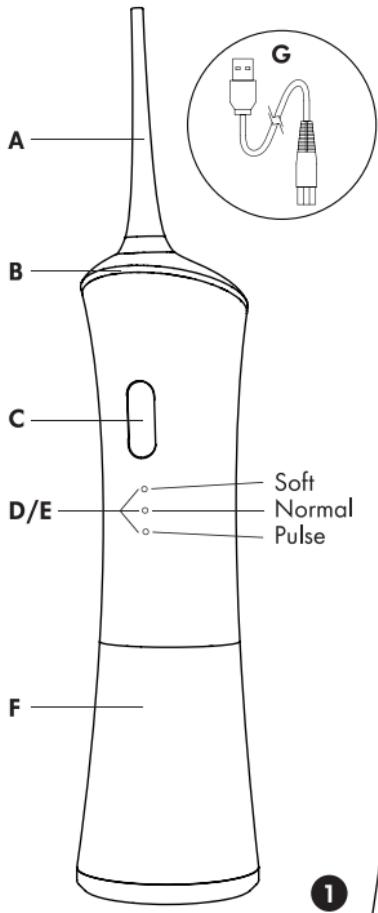PL: Tel.: +48717166648, www.mojadedra.pl

Product configuration

1x main body, 3x Nozzles, 1x USB Charge Cable, 1x user manual

Safety instructions

⚠ Dangerous

1. Please read the manual carefully before use, lest the improper operation cause injury or damage.
2. Do not disassemble this product privately, lest the product damage or cause electric shock accident.
3. Do not use more than 40 degrees of water into the tank, so as not to burn
4. If the product falls into the water when charging, please unplug the Power plug first.
5. If the power cord or plug is damaged, please replace the power adapter then use.
6. Avoid high temperature place this product, do not put the product into the fire.


⚠ Warning

1. Do not use the product for other purposes other than oral cleaning.
2. Do not use the nozzle pressure on the gums or teeth, so as not to damage the gums.
3. Do not immerse the product in the water to avoid the product failure.
4. Do not use this product in high temperature and humid environment.
5. Do not continue to use the power cord if the power cord is damaged.


⚠ Note

1. Do not use this product without installing a nozzle.
2. Some mouthwash may have a certain corrosive, do not put into the tank.
3. If the child or the operator does not want to use this product, to ensure that there is some supervision or guidance on the use.
4. Please keep children and babies out of reach, do not use as toys.
5. Patients with periodontal disease or in the last 2 months of oral surgery or suffering from oral problems, please consult the doctor before using.

Accessories:

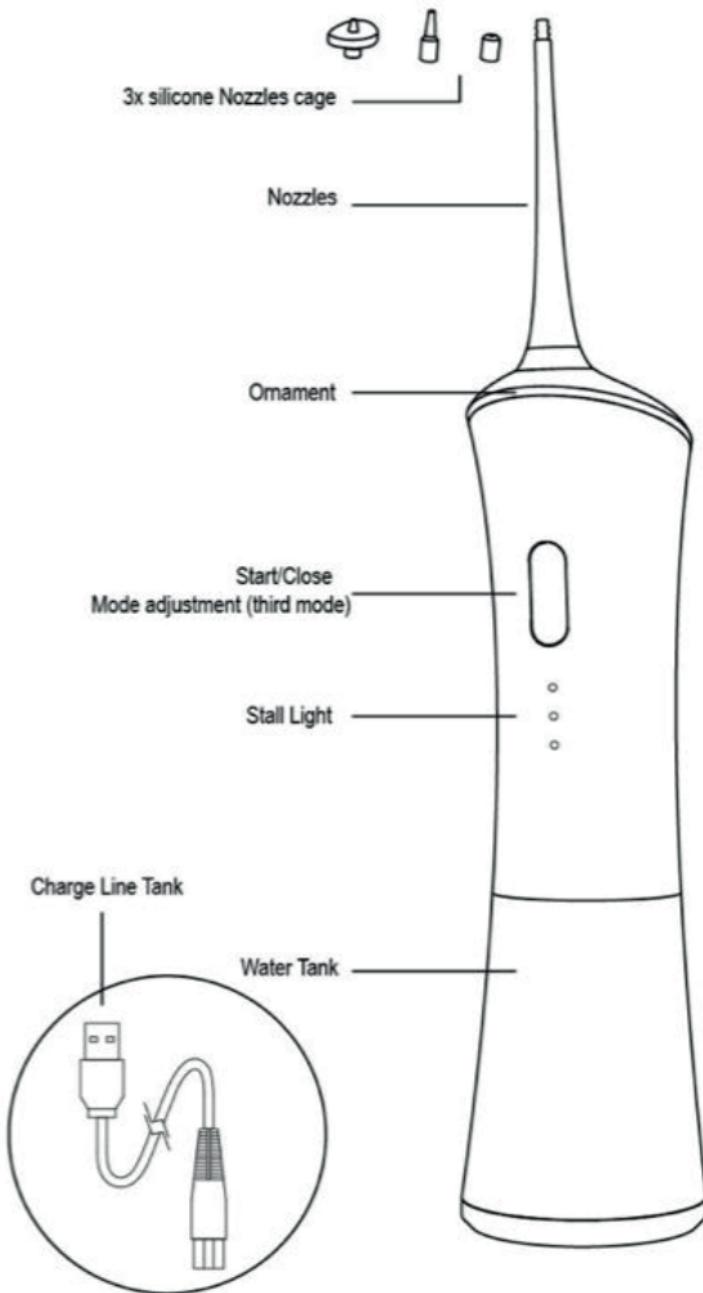Structure of the tooth-punching device

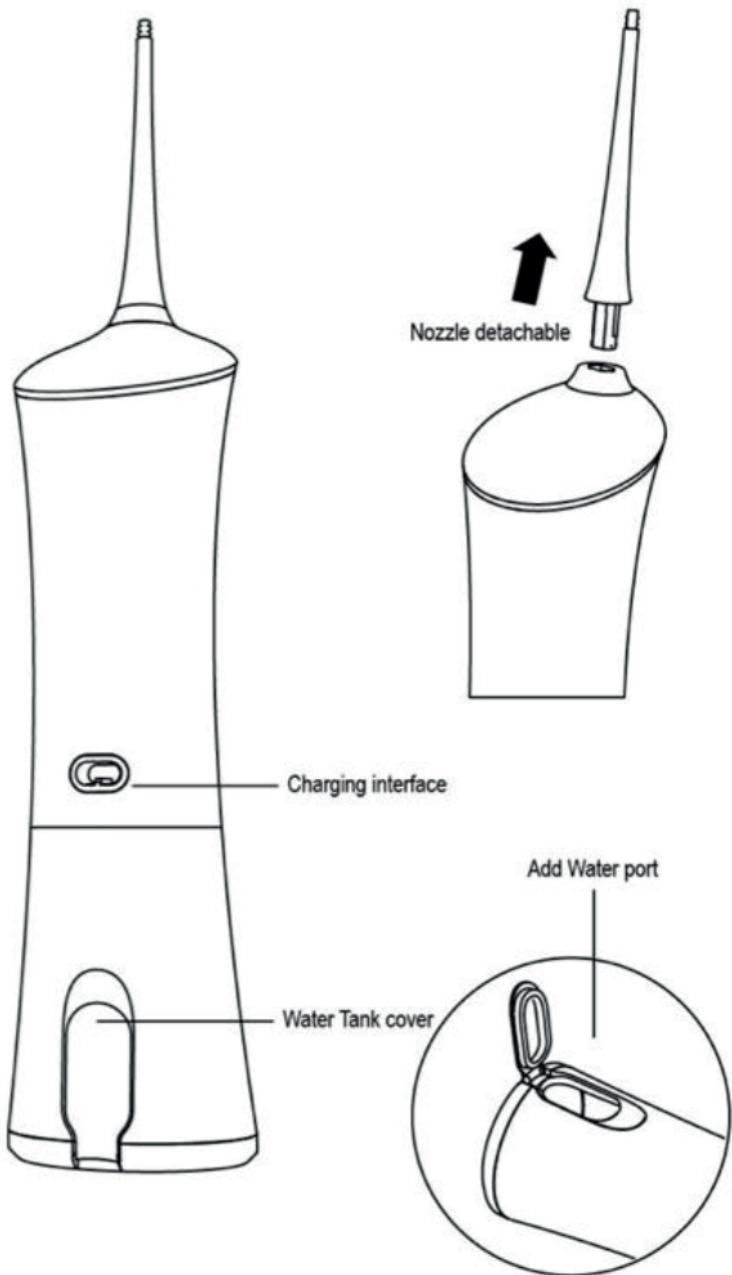Nozzle Silicone Sleeve

1. Short flat head: Suitable for most of the teeth, the water is moderate in thickness.
2. Long tip: Suitable for nozzle head depth on both sides of the teeth, but also suitable for exploring some of the more potholes.
3. crown Fitting head: Suitable for the head with teeth or gum surface, while moving edge flushing, and the flow is not easy to rebound outside.


Note:

Three kinds of nozzle silicone sleeve is to provide a more personalized use of experience, specific users can choose to use or not to use according to the actual situation.

Use steps

1. Add Water


Open Add water mouth clamshell, horizontal set this product, can directly put under the faucet to receive water. Or use a glass to put water, pour into the tank.

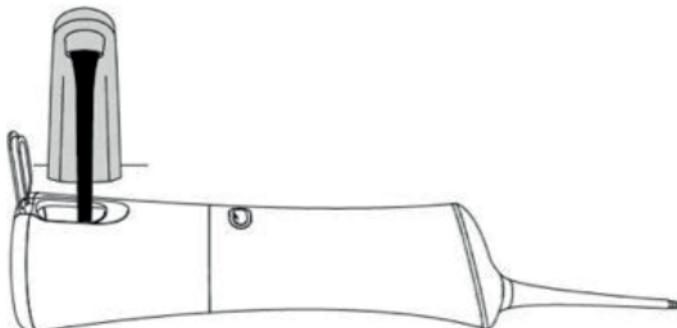2. Clean Teeth

Can choose to suit their own stalls, the proposal from the most low-grade start trial, and gradually find the most suitable for their own stalls.

When cleaning the teeth, extend the nozzle into the mouth, aim at the teeth or gums, and open the mouth slightly so that the flow of water flows smoothly.

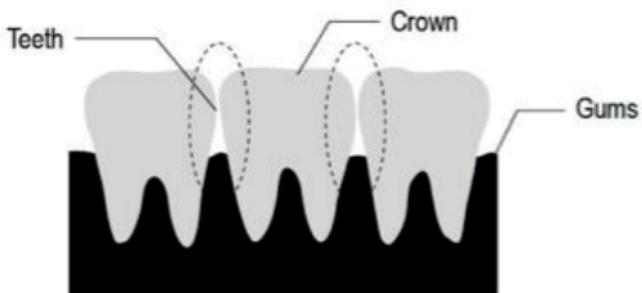⚠ Note

Do not use a strong file to continuously scour the gum damage, or dental bleeding place.

Please adjust the pulse flow intensity according to the condition of the oral cavity.

3. Discharge

Open Add water mouth clamshell, horizontal set this product, the remaining water dumped out.

4. Cleaning

⚠ Note

Do not put the whole machine under the faucet to clean or immersed in water cleaning, will lead to product failure or complete damage.

Please use cloth or other dry paper towels to wipe the teeth.

If not used for a long time, please make sure the water tank is dry and then stored.

Product Utility

- Strong pulsed currents can completely clean the mouth of bacteria, dental stains, tartar and food residues .
- Ultra-fine flow can clean gum groove and deep teeth without barrier, improve gum health .
- Effective removal of most stubborn dental plaque on the surface of teeth, the effect is much better than traditional dental floss .
- Quick cleaning of the mouth, whether at home or out of the dining, is an effective way to quickly restore clean mouth .
- For the correction of teeth, dental implants, dental veneer around the cleaning has a better effect . Help to develop healthy daily teeth cleaning habits and keep teeth clean at all times

Daily Maintenance

- If you need to clean this product, please use dry cloth or wet cloth to scrub the body, nozzle.
- Do not immerse the body in the water.
- Do not use more than 50 degrees of hot water cleaning.
- Do not bend, stretch and twist the nozzle hose.
- It is recommended that nozzles be replaced every 6 months For a long time, please clear the water in the tank and place it in the shady dry place.

Recovery

When discard the product, please do not throw it with the general living rubbish in the instrument, should Hand it over to a designated recycling center, which is beneficial to the environment.

FAQ

Q: Is it normal to use products that cause gums to bleed?

A: Gum bleeding is caused by inflammation caused by bacterial propagation between the teeth. Sticking to the teeth each morning and night can help most users to slow down inflammation and bleeding. If the regular use of the teeth after a few days, still bleeding condition, please consult your dentist

Q: How to judge the effect of the tooth punch?

A: The role of the punch is to remove dental plaque between the teeth, helping to improve the health of the gums. Dental plaque is not visible to the naked eye, so if the following symptoms appear, that means that the product has been effective: over time, the gingival bleeding situation is decreasing, the gingival tissue appears pink instead of red, the gums between the teeth appear to be sharp, rather than an inflamed arc-shaped. Usually after 2-4 weeks of use of the teeth, it will significantly improve the health of the gums. You can consult your dentist about your gum improvement next time you look at your teeth.

Q: Can the teeth be cleared of all the things between the teeth?

A: As part of the daily oral care program, the effect of the punch is to clear the plaque between the teeth, which should remove most of the food residue. However, it is not guaranteed to be used to remove large, and sticky food residues.

For more information, other FAQs and guidance videos, please go to official flagship store.

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