SYSTEMAT

BLACK HEADS FACIAL VACUMIZER **EL25931**

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1.Product Principle

EN

The pressure difference between the skin and the suction head is caused by the pneumatic pump in the product, sucked the blackheads, whiteheads and oil from the skin, so as to achieve the effect of cosmetic.

2.Accessory

Independent disposable suction head
No need to assemble, no need to clean, thrown away after use

big suction head

Larger suction, Suitable for large areas of forehead, Nosewing

flat suction head


sucked and quickly pulled out, repeat it one more time, with the effect of pulling compact skin.

small suction head

Suction is relatively weak, Suitable for thin, tender, fragile areas of skin

3.Product Parts

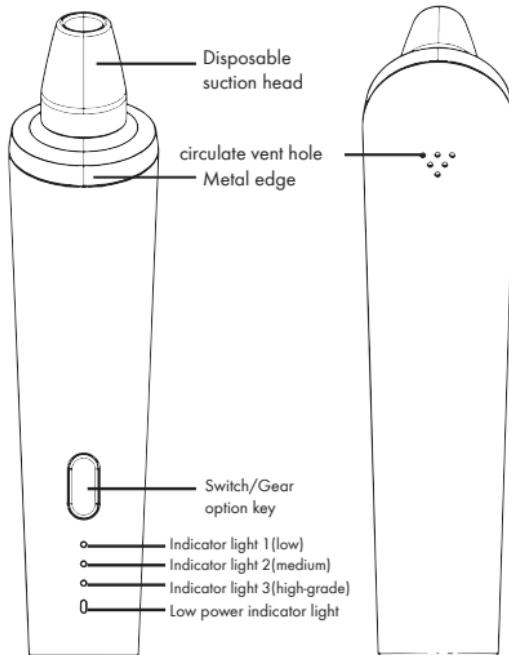ATTENTION

- 1) Far away from the water and wet place when charging.
- 2) When the power is less than 20%, the low light will flash, the low light bright when Charging, Position indicator light off when charging full.
- 3) 40minutes charging time,it can provide 50minutes use time when charging full.
- 4) Please charge between 0-40 degrees, avoid direct sunlight, avoid heat place when charging.
- 5) Please do not charging by other product usb data cable, so as not to damage the machine function.
- 6) Please provide the suitable power adapter when charging

4.Usage method

- 1)Please operate in the hand skin when you got this machine, Use on facial skin after mastering skills
- 2)Please do not use scrub or other horny skin care products before and after use this machine.
- 3)In order to avoid damage to the machine and your body, please remove all jewelry and contact lenses before use.
- 4)use it once per day(2minutes per once),excessive use may cause skin problems
- 5)Particular attention:do not stay in the same skin for a long time,so as not to cause swelling and discomfort
- 6)After nursing,skin appear red is normal phenomenon,it will become normal after 30 minutes later.
- 7)Recommended dry skin use it once a week,oily and other skin use it twice a week,2minutes per once,specific situation depends on individual skin conditions.

INAPPLICABLE SCOPE

- 1) Pimple acne has broken,scar on face
- 2) Eczema,abscess,acne,suppurative Diseases
- 3) Eyeball,eyelid and other delicate area of skin,so as not to damage the skin
- 4) The skin or the corresponding body parts has dermatitis or extreme sunburn
- 5) metal under the skin after plastic surgery
- 6) Rough skin

USE STEP

step 1: clean skin, hot compress

After clean skin. use facail steamer for 5minutes, or hot compress for 5 minutes, so that the pores open,black derived liquid can be match use .

step 2: select, install suction head

Choose your own need suction head, installed it on the machine

step 3: Boot Gear option

Press switch / Gear option to boot,according to you own skin condition select suction strength,

step 4 :Select area to start using

Select the area of the skin that you need to clean,slowly moving suction head on the skin, do not force the blackhead pressed on the skin or stay in the same area for too long time.

step 5: finish clean,shut down machine

The whole process takes 2 minutes. finish clean, press the switch / gear option for 2 seconds,shut down machine.

step 6 :cleansing skin, moisturizing skin

recommended finish clean, clean the face with cold water, Shrink pores and keep the face moist.

5.Matters needing attention

Please read carefully before use,do not follow the use of methods, product overheating, battery rupture, product damage, may bring harm to the human body.

- 1) the machine can not be immersed in water and washed.
- 2) dry your hand before use, so as water or other liquid far away from the back of vent hole,damaged the machine.
- 3) be sure to turn off the power after use the machine.
- 4) it will automatically shut down when the power is too low.
- 5) do not use for any other purpose than facial beauty,please stop use it immediately when you feeling pain or discomfort.
- 6) if the suction head is damaged, do not continue to use.
- 7) first-time use or skin sensitive should be shortened use of this product.

- 8) lead abnormal situation,immediately terminate use, and consult doctor in time.
9) prohibit children use,put it out of the reach of children.

6.Machine maintenance

- 1)do not put this product in a wet place.
- 2)do not let water to flow into the machine during clean.
- 3)do not use abrasive material to clean this product.
- 4)do not place heavy objects on the power cord and the product.
- 5)please shut down after use,in order to extend the service life of the battery,please full charge first,then storage this product for a month.

Special attention:

Do not disassemble, modify, repair the machine by yourself, please contact the after-sale maintenance while fault occurs.

7.Friendly reminder

Disposable suction head, No need to assemble, no need to clean, thrown away after use, you can buy it according to individual needs.

8.Fault and solution

Fault	Solution
failure to open the machine	<p>1.Check whether the power is insufficient, if the electricity is insufficient, please recharge</p>
succed to open the machine but suction strength is not normal	<p>1.Check whether Gear option, if not, select gear low medium or high-grade</p> <p>2.if you have chosen the gear, but the suction is not enough, please check whether the suction head is installed closely, whether leak or not,adjust and install suction head again.</p> <p>3.If the suction is too small, please check the LOW light, whether electricity shortage or not,please use after recharge</p> <p>4.if you have chosen the gear,suction head installation and contact are good,but machine does not work, please contact after-sale service department</p> <p>5.If machine is normal,but suction is not normal,for the sake of safety,please contact after-sale service department.</p>

Other unknown cause of failure,for the sake of safety,please contact after-sale service department.

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